

# THERMOSTIM



## 3 Treatment Modalities. 1 Tool.

Ready to deliver heat or cold therapy in less than 60 seconds.

**Ergonomic Design:** Specially designed grooves allow a close grip to control soft-tissue work.


**Soft-Tissue Mobilization:** The multi-surface head with flat and rounded surfaces, edges, and corners offers several options for soft-tissue mobilization.

### Heat/Cold & Electrotherapy:

Equipped with our patented technology, ThermoStim delivers heat or cold therapy from 35°F to 112°F in combination with any of 6 stim waveforms.



The **Dynatron ThermoStim Probe**, an accessory to the Dynatron Solaris® Plus, delivers 3 treatment modalities in 1 tool. Using ThermoStim, a practitioner can simultaneously deliver soft-tissue mobilization, heat or cold therapy, and electrical stimulation.

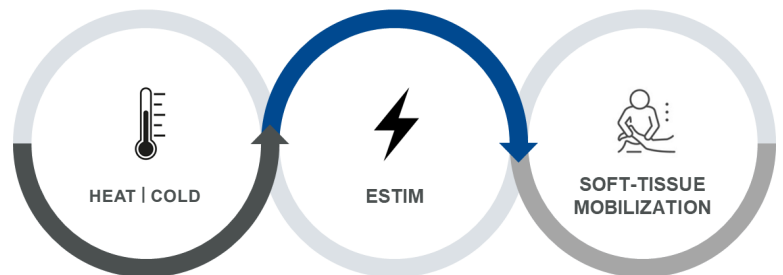

# THERMOSTIM



## ThermoStim in three easy steps.

**SET UP**  
Plug in cables, power on, and prep patient.

**1**

**2**

**3**

**SELECT**  
Select treatment parameters.

**START**  
Press start! ThermoStim will reach your selected temperature in <60 seconds.

### FUNCTIONALITY

The ThermoStim probe is designed to operate in conjunction with six separate stim modalities: Premod, Interferential, Biphasic, Russian, Micro-current, and High Volt. During a ThermoStim treatment, the probe becomes a mobile electrode as part of a stim treatment and delivers heat or cold to the patient at the same time.

### COMPATIBILITY

The ThermoStim probe (Part #: DTSP1) is compatible with all Dynatron Solaris® Plus Models: D715T, D716T, D718T, D719T.